
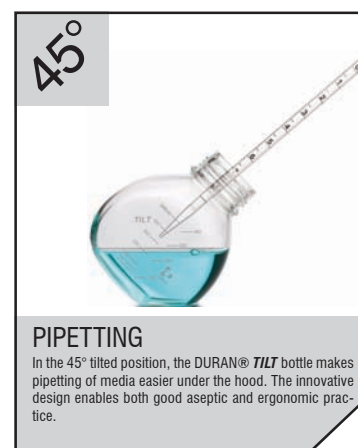


DURAN *TILT* TURNING USABILITY INTO SECURITY

Biosafety cabinets and clean hoods are essential for working with cell cultures under sterile conditions. Unfortunately, they make the handling of media bottles difficult and time-consuming for researchers. The DURAN *TILT* bottle changes everything. The *TILT* bottle has two positions: upright for filter sterilization or storage, and tilted at 45° for pipetting. Take a look at cell culture media from a new angle and discover a bottle system that eases handling and turns usability into security.

- Borosilicate 3.3 glass for reusability and lower whole life costs
- Space saving shape for fridge storage and water baths
- Low bottle height for easy handling in safety cabinets
- GL 56 mouth and screw cap - bigger opening for easy pipetting and safer handling
- 45° *TILT* position for pipetting without holding the bottle
- Tilted volume indicators for convenience
- Vacuum resistant bottle for safe filter sterilisation

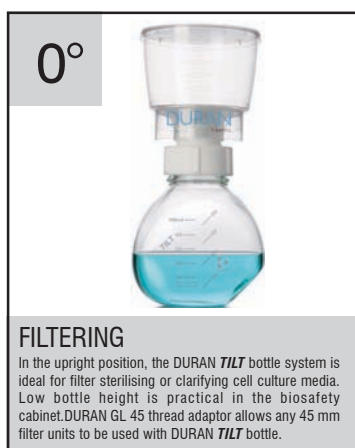

DURAN
Bottle top
vacuum
filtration
unit



DURAN
GL45
Thread
Adapter



DURAN
TILT
Protective
Light
Shield

DURAN *TILT* – MEDIA BOTTLE

Cat. No.	Description	Glass thread	Height (incl. cap) mm	Qty per Pack Pcs.
21891 4459	DURAN® <i>TILT</i> Bottle, clear glass Complete with white polypropylene screw cap 500 ml (supplied non-sterile)	GL 56	151	4

DURAN *TILT* – ACCESSORIES AND REPLACEMENT PARTS

Cat. No.	Description	Material Colour	Qty per Pack Pcs.
29243 5601	DURAN <i>TILT</i> Light Shield includes four GL 56 Bottle Tags (Orange, Yellow, Blue, Purple)	Silicone / White and mixed colours	4
29243 5626	DURAN GL 56 Bottle Tag Orange	Silicone / Orange	20
29243 5634	DURAN GL 56 Bottle Tag Yellow	Silicone / Yellow	20
29243 5659	DURAN GL 56 Bottle Tag Blue	Silicone / Blue	20
29243 5667	DURAN GL 56 Bottle Tag Purple	Silicone / Purple	20
29119 5601	DURAN GL 45 thread adaptor for bottle top vacuum filtration unit (GL 45 external, GL 56 internal)	PTFE / White	1
29401 5604	DURAN <i>TILT</i> GL 56 self-adhesive identification label for screw cap, 40 mm diameter	Polyester / White	1
29229 5602	DURAN <i>TILT</i> GL 56 screw cap	Polypropylene / White	10

DURAN BOTTLE TOP VACUUM FILTRATION UNITS

Cat. No.	Description	Membrane Pore Size μ mPES	Filter Dia mm	Thread mm	Qty per Pack Pcs.
29270 2818	DURAN 500 ml Funnel Only Gamma sterilized	0.10	90	45	12
29270 2826	DURAN 500 ml Funnel Only Gamma sterilized	0.20	90	45	12
29270 2842	DURAN 500 ml Funnel Only Gamma sterilized	0.45	90	45	12